

Model JSMI-04T 125 L Four Chamber Incubator



JSMI-04T 125 Liter Four Chamber Multi-Room Incubator provides two chambers in a compact unit share less footprint. Temperature of each chamber can be controlled independently. Temperature range from ambient +5 °C to +70 °C. Suitable for incubation and various test for food and packaging test on QA / QC labs.

MAIN FEATURES

PERFORMANCE

- Control Temperature Range from ambient +5 °C to +70 °C
- Forced air circulation for excellent thermal accuracy and uniformity
- Aero-dynamic chamber designed with circulation blower to achieve superior thermal uniformity
- Separate heating compartment located on back wall for contamination free and maintenance free

CONSTRUCTION


- Durable epoxy powder coated metallic casing
- Corrosion resistant Stainless Steel 304 chamber with seamless round cornered edge allows easy cleaning and prevent leakage
- Tempered safety inner glass door for easy observation of samples

CONTROLLER




-3 ZONE PID ALGORITHM

Usable temperature range divided into three zones (LOW/MIDDLE/HIGH) and each ZONE has its individual PID parameters optimized to the range, and it provides precise and superior temperature control characteristic to whole temperature range.

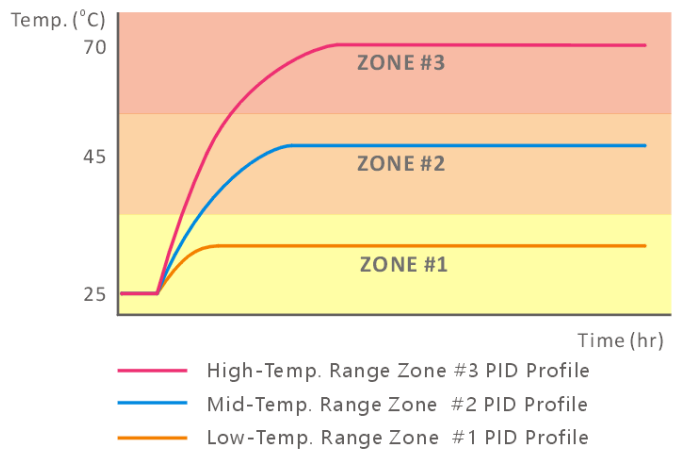


-3 POINT CALIBRATION FUNCTION

A perfect calibration tool to calibrate chamber at THREE DIFFERENT TEMPERATURE POINTS to adjust multiple temperature bias to increase the temperature accuracy of WHOLE TEMPERATURE RANGE.

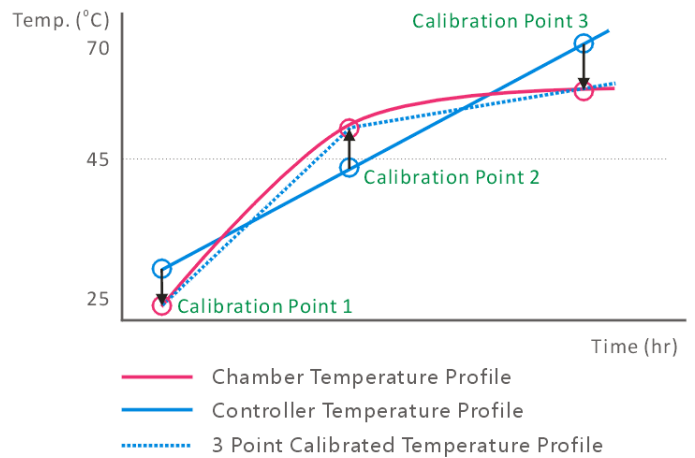

- 1-POINT OFFSET FUNCTION

Provides simple and easy calibration tool to calibrate chamber at ONE SPECIFIC TEMPERATURE POINT in frequent use.

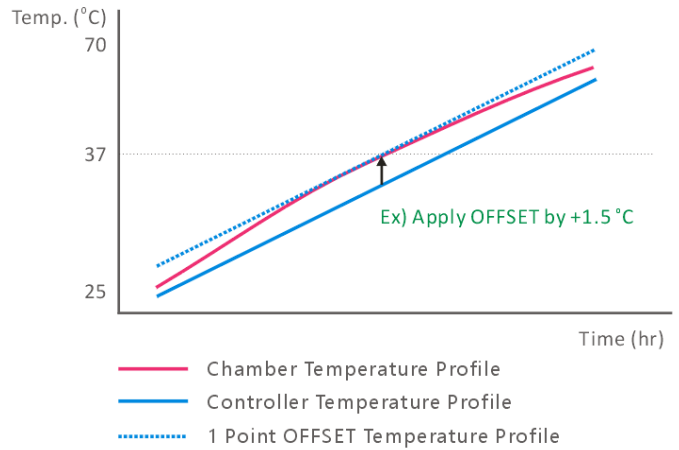

CONTROLLER

-DIGITAL MICROPROCESSOR PID TEMPERATURE CONTROLLER

Integrated Digital Microprocessor Control System

-GRAPHIC LCD DISPLAY

Back Light Graphic LCD Display provides user-friendly operating information

-UV COATED POLYCARBONATE TOUCH PANEL

Equipped with Capacitive LED Lighting Touch Button provides sensitive touching operation and durability

-ENCODER TYPE JOG DIAL

Push to Start and Stop, Rotary turn to change values provides easy and simple user interface

-0.1 °C CONTROL RESOLUTION

-Auto-Stop Timer 1 min ~ 99 hour 59 min and alarm

-Auto-Start Timer 1 min ~ 99 hour 59 min and alarm

-Temperature and Timer Memory Function in frequent use

-3-ZONE PID ALGORITHM

-3-POINT CALIBRATION

-1-POINT OFFSET CALIBRATION

-POWER and TIMER RESTORE FUNCTION in case electric power outage

-CLASS A PT-100 Ω Sensor

-Delay starter to protect compressor.

-RS-232C PC COMMUNICATION INTERFACE [OPTION]

CONVENIENCE

-Outer door with self-adhesive magnetic packing

-Inner glass door with high density silicone gasket and air-tight door lock

-Removable and height adjustable shelf by 25mm spacing

SAFETY

- PT SENSOR SHORT OR DISCONNECTION AUDIO VISUAL ALARM
- DUAL OVER TEMPERATURE CUT-OFF BY CONTROLLER AND INDEPENDENT THERMOSTAT
 Controller cut-off heater, compressor and AUDIO VISUAL ALARMS
 ALARM HIGH : +5 °C above set temp.
 ALARM LOW : -5 °C below set temp.
 Independent analog thermostat cut-off heater above set temperature
- DOOR OPEN AUDIO VISUAL ALARM
- OVER CURRENT CUT-OFF : Electric Circuit Breaker

SPECIFICATIONS

JSMI-04T

Dimension Data

Chamber Volume	4 chamber x 125 L
Inner Dims. (W x D x H mm)	500 x 500 x 500
Outer Dims. (W x D x H mm)	1365 x 830 x 1815
No. of Shelf Positions	16 ea
No. of Shelf included	8 ea

Performance Data

Heating System	Forced Air Convection
Cooling System	none
Temp. Range	amb.+5 °C ~ +70 °C
Accuracy at +20.0 °C	±0.2 °C at +37.0 °C
Uniformity at +20.0 °C	±1.0 °C at +37.0 °C

Controller

Processor	Digital Microprocessor PID Temperature Controller
Display	Back Light Graphic LCD Display
Panel	UV Coated Polycarbonate Solid Touch Panel
Button	Capacitive touch LED Lighting Touch Button
Interface	Push and Rotary encoder type Jog Dial
Control Logic	3 Zone PID Control
Calibration	1 Point Temperature OFFSET 3 Point Temperature Calibration
Resolution	0.1 °C
Auto-Stop Timer	1 min. ~ 99 hr 59 min
Auto-Start Time	1 min. ~ 99 hr 59 min

Technical Data Model JSMI-04T

Sensor	CLASS A PT-100 Ω Sensor
Power Restoration	Power and Timer Restoration Function in case of power outage
PC Communication Interface	RS-232C PC COMMUNICATION INTERFACE [OPTION]

Safety Data

Safety Device	Sensor short or disconnection AUDIO VISUAL ALARM Controller cut-off heater, compressor and AUDIO VISUAL ALARM HI-ALARM : +5 °C above set temp. LO-ALARM : -5 °C below set temp. Independent analog thermostat cut-off heater above set temperature Door open AUDIO VISUAL ALARM
---------------	--

Material Data

Body	Epoxy Powder Coated Steel
Chamber	Stainless Steel 304 (AISI304)
Inner Glass Door	Tempered Safety Glass, 5 mm thick

Electrical Data

Rated voltage	220 \pm 10% VAC
Power frequency	50 / 60Hz
Power Cord	Two for 4 chambers
Heater output	4 x 200 W
Nominal power total	880 W 4.0 amp
Nominal power per plug	440 W 2.0 amp
Phase	1 Phase
Plug configuration	CEE 7/7 Schuko or BS 1363

Usable Condition

Temp. & Humidity	5°C - 40°C, < 80% RH
------------------	----------------------

Packaging Data

Type of packaging	Wooden Crating
Shipping Dims. (WxDxH mm)	1510 x 960 x 2000
Net Weight	220.0 kg
Shipping Weight	350.0 kg

ORDERING INFORMATION

Art. No.	Descriptions	Voltage
JSMI-04T	125 L Four Chamber Incubator	220 VAC 1~ ph 50/60Hz

OPTIONAL ACCESSORIES

Art. No.	Descriptions
PHS-A06	Perforated Shelf, Stainless Steel 304, W470 x D480 mm for JSMI-04
PHS-C05	PVC Coated Wire Shelf, W470 x D480 mm for JSMI-04
PHS-A99	Shelf Support, Set of 4 ea

JS RESEARCH INC.

Gongju-City, Korea

TEL +82 41 856 3677

FAX +82 41 856 3679

info@jsresearch.co.kr

www.jsr.kr